
## **Open-Source Syringe Pump**

(developed by Andrey Samokhin in 2018)

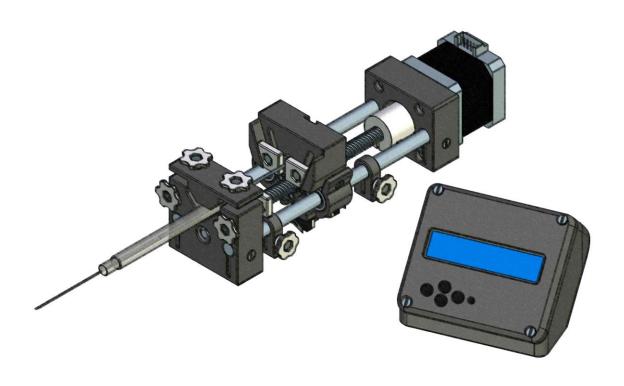

## **Assembly Instructions**

(Version 1.0a)

STL-files and "firmware" are available at <a href="http://www.mass-spec.ru/projects/diy/syringe\_pump/eng/">http://www.mass-spec.ru/projects/diy/syringe\_pump/eng/</a> Project files are distributed under CC BY-SA 3.0 license



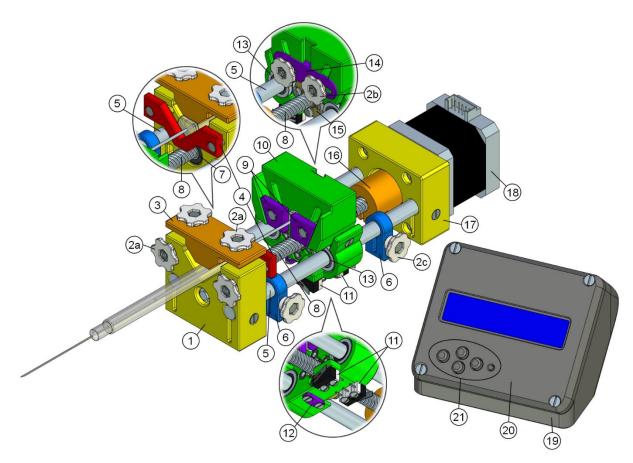


Fig. 1. General view.

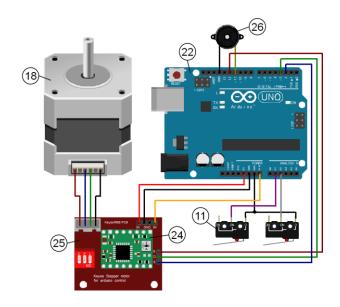


Fig. 2. Electrical circuit.

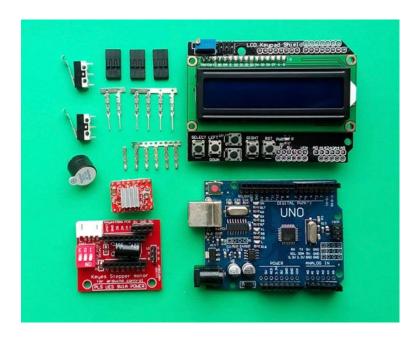


Fig. 3. Electronic components.

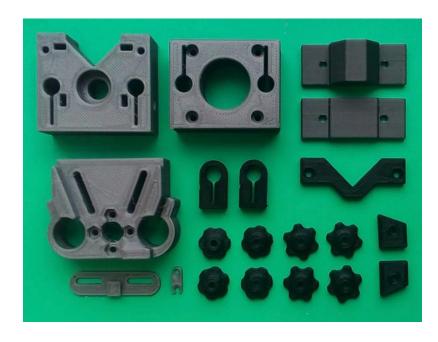


Fig. 4. Plastic parts made on 3D printer.



Fig. 5. Plastic case made on 3D printer.




Fig. 6. Guide rods, lead screw with nut, metal rigid coupling, bearings, stepper motor, screws and nuts.

Table 1. Part list.

PART#	PART NAME	QTY
1	Front Support	1
2a	Hand Knob (9 mm)	4
2b	Hand Knob (8 mm)	2
2c	Hand Knob (7 mm)	2
3	Top Syringe Holder	1
4	Side Syringe Holder	1
5	Guide Rod (D=8 mm)	2
6	Limit Stop	2
7	Bearing 688ZZ	1
8	Lead Screw (T8)	1
9	Plunger Holders	1
10	Carriage	1
11	Endstop Switch	2
12	Wire Holder	1
13	Linear Bearing (LM8SUU)	2
14	Slider	1
15	Lead Screw Nut (T8)	1
16	Coupling	1
17	Back Support	1
18	Stepper Motor (NEMA 17)	1
19	Case (Base)	1
20	Case (Cover)	1
21	Buttons	1
22	Arduino UNO	1
23	LCD Keypad Shield	1
24	Stepper Motor Driver (A4988)	1
25	Control Shield for A4988	1
26	Active Buzzer (D=12 mm)	1
27	Nut M3	20
28	Countersunk Screw M3x12	4
29	Countersunk Screw M3x14	4
30	Countersunk Screw M3x20	4
31	Pan Head Screw M3x16	4
32	Bolt M3x20	2
33	Bolt M3x25	2
34	Bolt M3x35	4
35	Self-tapping Screw 2x8	6
36	Self-tapping Screw 2x10	8
37	Dupont Connector (Female)	6
38	Dupont Connector (Male)	6
39	Dupont Housing (3P)	3
40	Connector Wires	_
41	Power Supply (12 V, 2 A)	1

Step 1.1 (assemble according to the drawing)


<u></u>		
PART #	PART NAME	QTY
8	Lead Screw (T8)	1
16	Coupling	1
18	Stepper Motor (NEMA 17)	1
27	Nut M3	2
28	Countersunk Screw M3x12	2

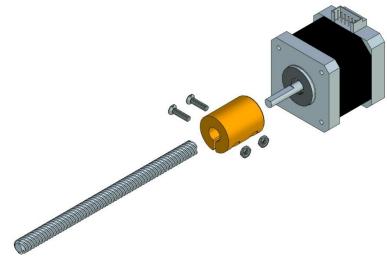

<sup>\*</sup> metal rigid coupling can be used instead of 3D printed part

Step 1.2 (assemble according to the drawing)


PART#	PART NAME	QTY
17	Back Support	1
27	Nut M3	2
30	Countersunk Screw M3x20	4

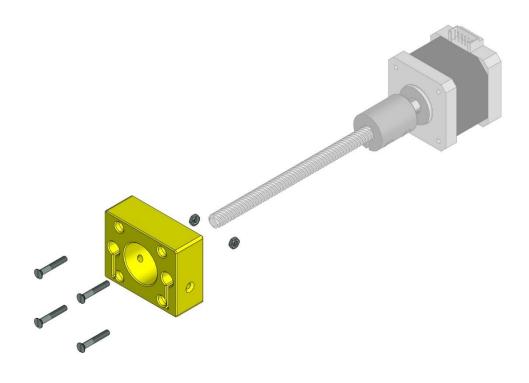

<sup>\*\*</sup> flexible beam coupling is not recommended to use


Step 1.3 (assemble according to the drawing)

	_ `	
PART#	PART NAME	QTY
5	Guide Rod (D=8mm)	2
29	Countersunk Screw M3x14	2

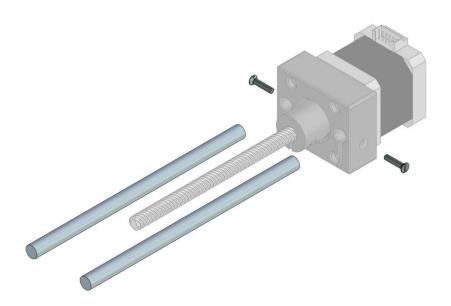

**Step 1.4** (assemble according to the drawing)


otop :::: (accernate acceraing to the ar		
PART#	PART NAME	QTY
10	Carriage	1
11	Endstop Switch	2
12	Wire Holder	1
13	Linear Bearing (LM8SUU)	2
15	Lead Screw Nut (T8)	1
27	Nut M3	6
28	Countersunk Screw M3x12	2
31	Pan Head Screw M3x16	4
35	Self-tapping Screw 2x8	2
36	Self-tapping Screw 2x10	4

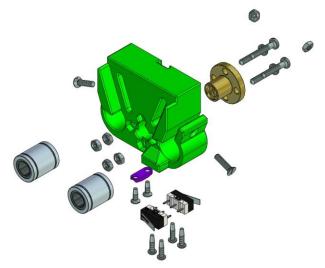

<sup>\*</sup> before installing the endstop switches solder the wires (refer to step 2.3)


Step 1.5 (assemble according to the drawing)

- 3		_ `	
	PART#	PART NAME	QTY
	2b	Hand Knob (8mm)	2
	9	Plunger Holders	1
	14	Slider	1
	27	Nut M3	2
	34	Bolt M3x35	2

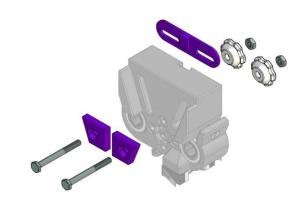


Step 1.6 (assemble according to the drawing)

PART#	PART NAME	QTY
2c	Hand Knob (7mm)	2
6	Limit Stop	2
27	Nut M3	2
32	Bolt M3x20	2

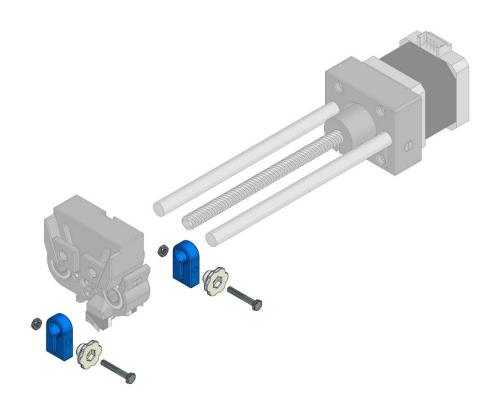


Step 1.7 (assemble according to the drawing)

PART#	PART NAME	QTY
1	Front Support	1
2a	Hand Knob (9mm)	4
3	Top Syringe Holder	1
4	Side Syringe Holder	1
27	Nut M3	4
33	Bolt M3x25	2
34	Bolt M3x35	2

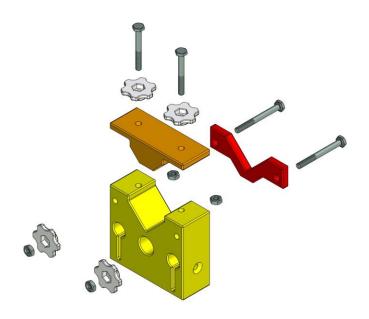

Step 1.8 (assemble according to the drawing)

	_ \	
PART#	PART NAME	QTY
7	Bearing 688ZZ	1
27	Nut M3	2
29	Countersunk Screw M3x14	2

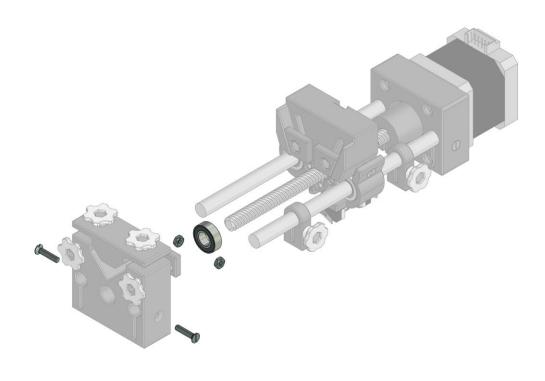


## Step 2.1 (if necessary, change the basic settings in the Arduino sketch)

## Step 2.2 (upload sketch to Arduino)

Detailed instructions are available at this link.

Step 2.3 (solder wires and crimp connectors, refer to electrical circuit)

PART#	PART NAME	QTY
22	Arduino UNO	1
26	Active Buzzer (D=12mm)	1
37	Dupont Connector (Female)	6
38	Dupont Connector (Male)	3
39	Dupont Housing (3P)	2
40	Connector Wires	_

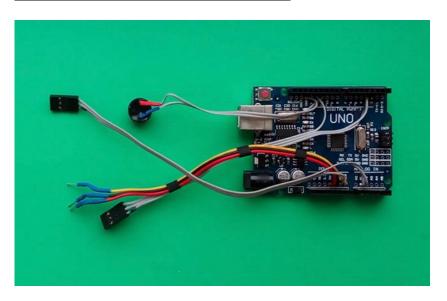

Step 2.4 (solder wires and crimp connectors, refer to electrical circuit)

PART#	PART NAME	QTY
11	Endstop Switch	2
38	Dupont Connector (Male)	3
39	Dupont Housing (3P)	1
40	Connector Wires	_




Step 2.5 (plug A4988 driver into the shield and set all DIP switches to ON)

PART#	PART NAME	QTY
24	Stepper Motor Driver (A4988)	1
25	Control Shield for A4988	1




<u>Step 2.6</u> (set output current limit on A4988) Detailed instructions are available at <u>this link</u>.

<u>Step 2.7</u> (connect all electronic components together; connect stepper motor and endstop switches; plug in 12 V power supply; make sure that pump works)

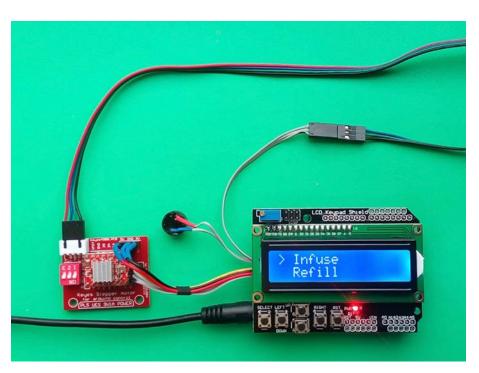


Step 2.8 (place active buzzer, Dupont housing, and shield with A4988 into the case)

PART#	PART NAME	QTY
19	Case (Base)	1
35	Self-tapping Screw 2x8	2

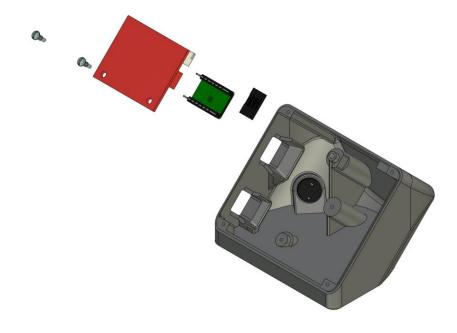


Step 2.9 (place Arduino UNO into the case)

PART #	PART NAME	QTY
35	Self-tapping Screw 2x8	2

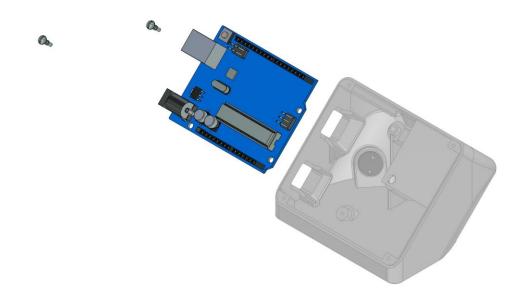

Step 2.10 (plug LCD Keypad Shield into Arduino UNO)

<u> </u>				
PART#	PART NAME	QTY		
23	LCD Keypad Shield	1		

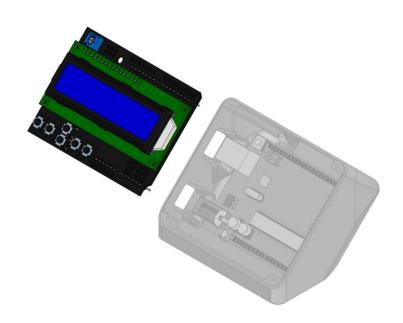

Step 2.11 (install buttons and top cover)

PART#	PART NAME	QTY
20	Case (Cover)	1
21	Buttons	1
36	Self-tapping Screw 2x10	4



